

Stools & Cha



Designed and manufactured to exceed your expectations...

It begins with a wide, five-leg caster base that provides support and stability. The star-shaped base allows closer accessibility to the patient while heavy duty, dual-wheel, hooded casters enable stools to maneuver quietly over all types of surfaces.

Adding to the strength and stability of our stools is the solid steel machine screw on manual models, or air cylinder on pneumatic models. The air cylinder further acts as a built-in shock absorber, taking the abuse of a hard day's work instead of you, and presents a simple and efficient method for adjusting the seat height; ideal for multi-physician practices.

Comfort is provided in the form of a wide, amply padded seat section constructed of dense foam. A beadless upholstery seam eliminates the seat edge from digging into your legs. Add a backrest to promote better posture and support to leave you feeling less fatigued at day's end.

With many choices, colors and accessories, there is certain to be a stool just right for you!

276 Air Lift Stool

Our best selling stool! The 276 has become the industry standard in healthcare facilities.



Classic Series Stools

The healthcare industry's most demanded stools, and the foundation of our seating product line. The Classic Series features a polished aluminum/black accented five-leg caster base and hood, and a hand ring or foot lever height control on pneumatic models. The seat height of the 274 (without back) and 275 (with back) adjusts manually while the 276 (without back) and 277 (with back) adjusts pneumatically. The 278 (foot lever) and 279 (hand ring) offer pneumatic height adjustment and a surgical armrest to stabilize the surgeon's arm during delicate procedures.



Stools and Chairs

S P E C I F I C A T I O N S

Classic Series Stools
 Seat Height:
 18" - 24.5" (274-278)
 19" - 25.5" (279)
 Seat Diameter: 16.25"
 Seat Cushion Thickness: 3.25"
 Back Cushion Dimensions:
 14" W x 7.5" H (275, 277)
 Caster Base Diameter: 21.5"
 Weight Capacity:
 350 lbs. (274, 276)
 225 lbs. (275, 277, 278, 279)

Value Series Stools

Seat Height:
 18" - 24.5" (270-273)
 Seat Diameter: 16.25"
 Seat Cushion Thickness: 3.25"
 Back Cushion Dimensions:
 14" W x 7.5" H (271, 273)
 Caster Base Diameter: 23"
 Weight Capacity:
 300 lbs. (270, 272)
 225 lbs. (271, 273)

Classic and Value Series Options

Locking casters
 Soft rubber casters
 Glides
 6" Lab height extension
 Foot ring

425/427 Air Lift Physician Stools

Seat Height:
 17.5" - 24" (425)
 17" - 23" (427)
 Seat Cushion Dimensions:
 16.5" W x 17" D
 Seat Cushion Thickness: 3"
 Back Cushion Dimensions:
 13" W x 9.5" H
 (with back kit accessory)
 Caster Base Diameter: 23"
 Weight Capacity: 225 lbs.

281 Blood Drawing Chair

Seat Height: 20"
 Seat Cushion Dimensions:
 20" W x 16.5" D
 Seat Cushion Thickness: 2"
 Back Cushion Dimensions:
 20" W x 16.5" H
 Inside Arm Width: 21"
 Outside Arm Width: 32"
 Storage Drawer Capacity:
 8.5" W x 15.25" D x 4.5" H
 Weight Capacity: 400 lbs.

680 Side Chair

Seat Height: 17"
 Seat Cushion Dimensions:
 18" W x 19.5" D
 Seat Cushion Thickness: 2"
 Back Cushion Dimensions:
 18" W x 13.5" H
 Inside Arm Width: 17.5"
 (on model with arms)
 Outside Arm Width: 23.5"
 (on model with arms)
 Weight Capacity: 400 lbs.


Upholstery Colors:




Ritter: Navy, Dusty Blue, Moss, Clay, Pebble Grey, Shadow, Perfect Plum and Black



Midmark Premium Solid: Blueberry, Silver Sage, Tea Green, Fossil Grey, Cashmere Blue, Terra Cotta and Black

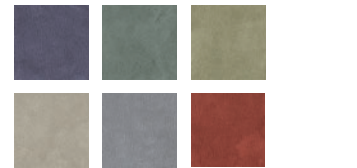

Midmark Premium Pattern: Belagio, Siena, Tuscany, Firenze, Capri and Treviso

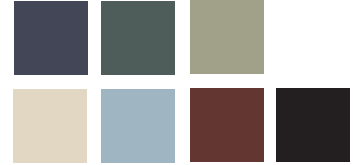

Midmark Ultra-Premium: Ultra-Pacific Blue, Ultra-Spruce, Ultra-Olivine, Ultra-Arctic, Ultra-Sky, Ultra-Canyon and Ultra-Raven Wing

*There will be an additional charge for Midmark Ultra-Premium.

Midmark is an ISO 13485 and ISO 9001 Certified Company.

For more information or a demonstration, contact your Ritter or Midmark Dealer, or contact Midmark Corporation, Versailles, Ohio. Call: 1-800-MIDMARK Fax: 1-800-365-8631 midmark.com

© Midmark Corporation 2005
 Products subject to improvement changes without notice.
 Litho in U.S.A. 007-0558-00 Rev. D1 (3/08)



Classic and Value Series Backrest



Classic Series Glides



Value Series Glides



Classic and Value Series Locking Casters



Classic Series Soft Rubber Casters



Value Series Soft Rubber Casters



Classic Series 277 shown with optional 6" lab height extension and foot ring.

Ritter®
 by MIDMARK