



500

Pediatric Crib

Stretcher



PEDIATRIC STRETCHER/CRIB

3177

Function as Appealing as its Form

For the care of your young patients, Pedigo's 500 Pediatric Crib Stretcher with its full function capabilities represents the standard of excellence in hospital cribs. The aesthetically pleasing design makes it perfect for any facility.

As a pediatric transport stretcher, the 500 rolls easily on 5" oversized casters. Two of the four casters lock to secure the stretcher while a third steering caster may be engaged for increased control and maneuverability.

Optional IV poles allow for bottles or bags to be suspended inside or outside the stretcher area.

The unique rounded corners reduce the overall size of the 500, yet provide the same amount of useful area on the 26½" x 55¾" mattress. An optional

utility shelf located at the base may be used for storing blankets or the child's personal items.

Total Accessibility

The two-post support system provides 97% access to the patient by the medical team.

Fixed Security Top

With the optional translucent security top in position, the 500 becomes a full function pediatric crib. The molded top has an exclusive quick disconnect for fast, easy removal.




Optional Departmental Color Coding

Identify your 500s with the departmental color of your choice on each end.

A True One-Handed Crib

Dual-action release handle is childproof and easy to operate, even with a child in one arm. Side rails have two-inch bar spacing and a special epoxy finish that is warmer to the touch than chrome.

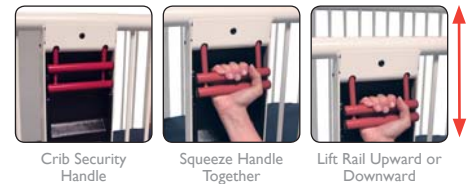

For complete containment of the child, simply raise the counter balanced side rails to enclose the crib. Because the rails are counter balanced, there will be no sudden drop when released, preventing injury to the child or provider.

500 STANDARD FEATURES

NEW!

- Two 5" locking casters and one steering caster for increased maneuverability
- Trendelenburg/Reverse Trendelenburg
- One handed rail operation
- White powder coat finish
- Unique rounded corners
- Rails with 2" bar spacing
- Two IV pole receptacles
- Full function fowler
- Litter Top: Solid pan
- Color decals: Blue

 **LATEX-FREE**



SPECIFICATIONS

Model Number	500 Stretcher Crib
Patient Surface	26½" x 55¾" (67.3 cm x 141.6 cm)
Overall Length	60" (152.4 cm)
Overall Width	30¼" (76.8 cm)
Height	
Fixed Litter Top	32¼" (81.9 cm)
Top of Crib	59½" (151.1 cm)
With Security Top	74" (188 cm)
Litter Positioning	
Fowler Angles	0°, 20°, 30°
Trendelenburg	0°, 5°, 10°
Reverse Trendelenburg	-5°
Side Rail Intermediate Stop	13¾ (39.9 cm)
Casters	5" (12.7 cm)
Warranty	One year limited

OPTIONS

5103001	Utility shelf	5114001	Oxygen tank holder
5104002	Mounting kit for traction adapter and croupette adapter	5116001	IV Pole
5106001	Infusion pump adapter	5860002	Mattress: Foam pad with Staph CHEK™ cover
5107001	Windowed Security top with 6 mounting clips (will accommodate a mist tent)	5117001	Security top, translucent, with 6 mounting clips (will accommodate a mist tent)
5108001	Croupette humidifier adapter (requires 5104002 for mounting)	5119001	Traction adapter, pair (includes mounting kit)
5109002	Crib bumper pad	5817002	Pressure care pad, 4"
		Color decals	